

LR Delivery Unit Installation



Parts Bag Contents

Stop Pins 2

Special Tools

Level

Torque Wrench

Attention!

Do not remove tape or packaging until Step 2 is complete.

CE 0088

Step 3: Install two stop pins.


Tip: Move stop pins under trim ring, then down into holes.

Step 1: Remove connection box cover.

Step 2: Attach arm to casting.
Do not tighten bolt!

Step 4: Screw bottle into cap.

DA134100i

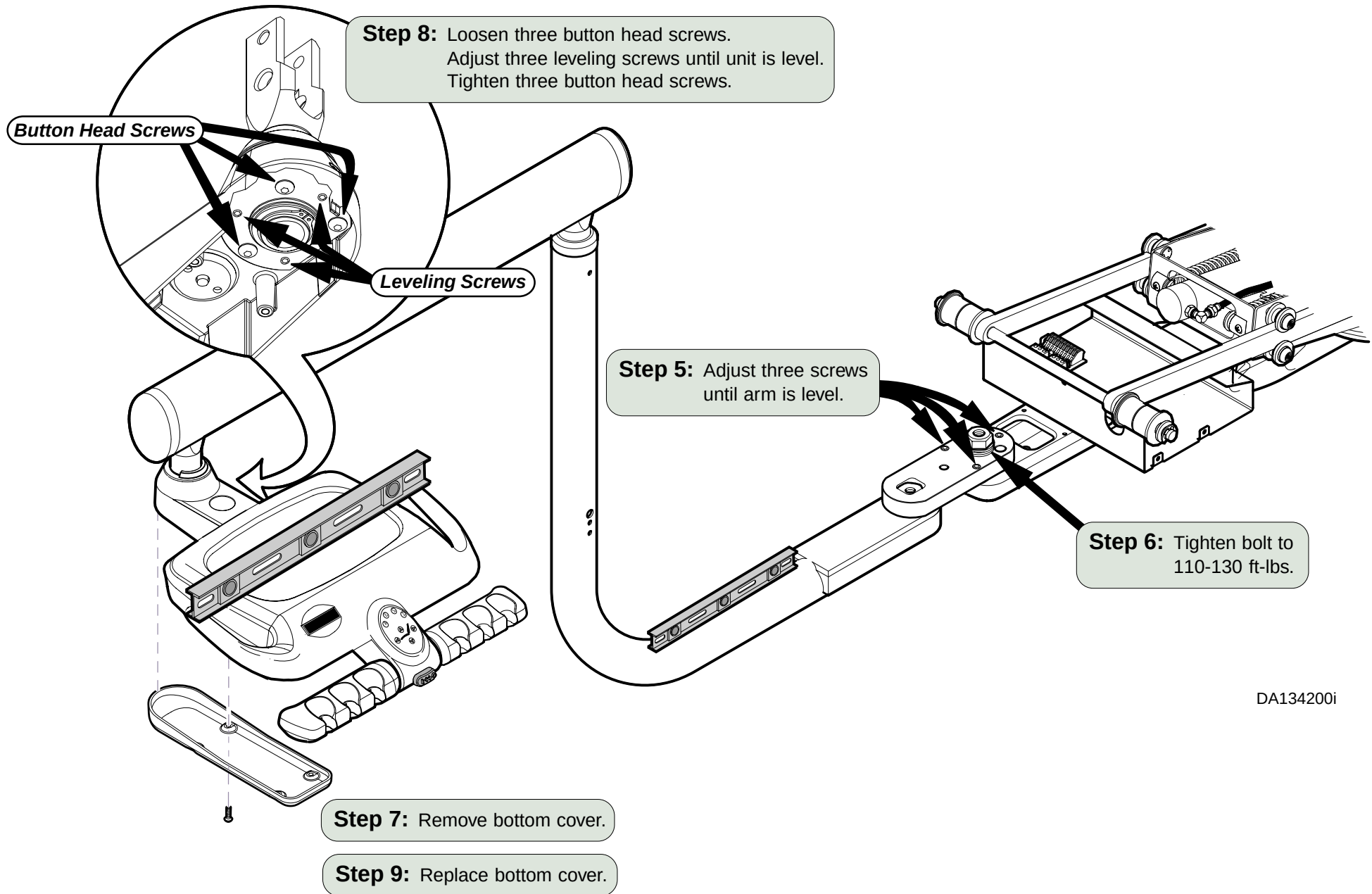

Attention!

Connect wire harnesses and tubing before placing umbilical in notch of connection box.


Step 10: Connect tubing as shown.

Accessories

Step 11: Connect accessory wire harness to J4 on PC Board.

Chair Position Touch Pad

Step 12: Connect touch pad wire harness to J2 on PC Board.

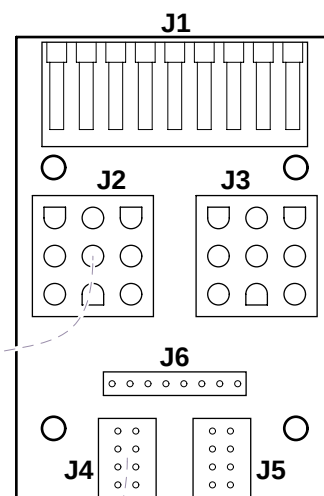

Step 14: Install connection box cover.

Step 13: Position umbilical in notch of connection box.

DA134300i

Attention!

Perform these steps after installation of entire operatory is complete.

